

# No Bubble TURP Irrigation set

## For Uninterrupted Resection and Precision Infusion

The No-Bubble TURP Irrigation set is the result of the input from a number of skilled experienced urology nurses and medical practitioners.

Made of extra soft material, the tubing rests snugly on the drapes and the drip chambers are effortlessly filled with irrigation fluid.

The spikes are positioned close to the bag for easy de-aeration, and both lines are clearly colour coded and easily separated, allowing for the pinch clamps to be moved just under the drip chambers. This makes it easily visible for when patients undergo continuous bladder irrigation in the post-operative period.

The colour coded pinch clamps are large enough to be opened and closed using only one hand and are placed in the sterile field during operation to allow the scrub nurse or surgeon to change from one irrigation bag to another without interrupting the procedure or the need for the circulating nurse's instant presence.

The drip chambers are equipped with water dispersers to reduce kinetic energy and number of bubbles during rapid irrigation.

The entire set is made of DEHP-free material to reduce the impact of disposal on the environment.



The Spike is attached directly to the drip chamber making piercing of the bag easier



Bubbles and fogging are eliminated by perforations in the inner dispersing tube that flow evenly along the walls, while maintaining a good flow rate




The softer drip chamber is easier to squeeze even if irrigating fluid is cold




Colour coded clamps help to identify the bag that is being drained

## Intra-Operative Non-Sterile

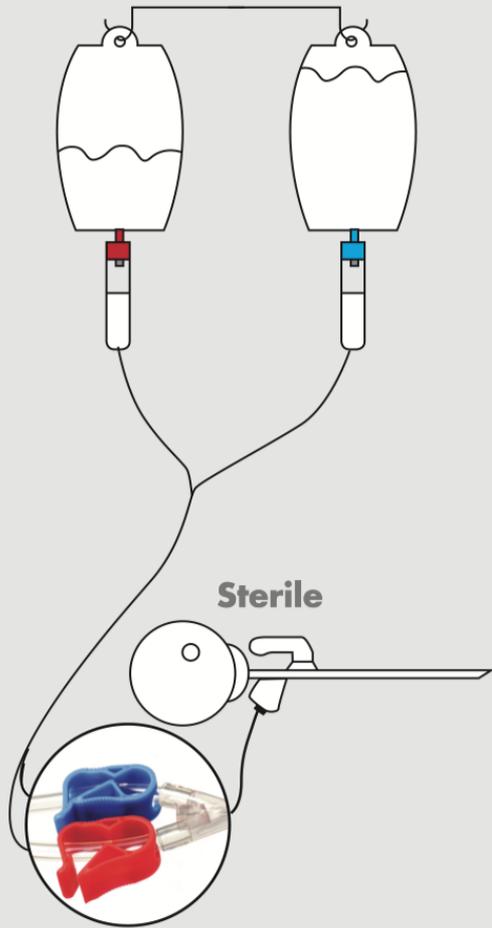

The clamps are located in the sterile field. The switch from an empty to a full bag is done either by a scrub sister or the surgeon. This allows uninterrupted resection and saves operating time.

## Post-Operative

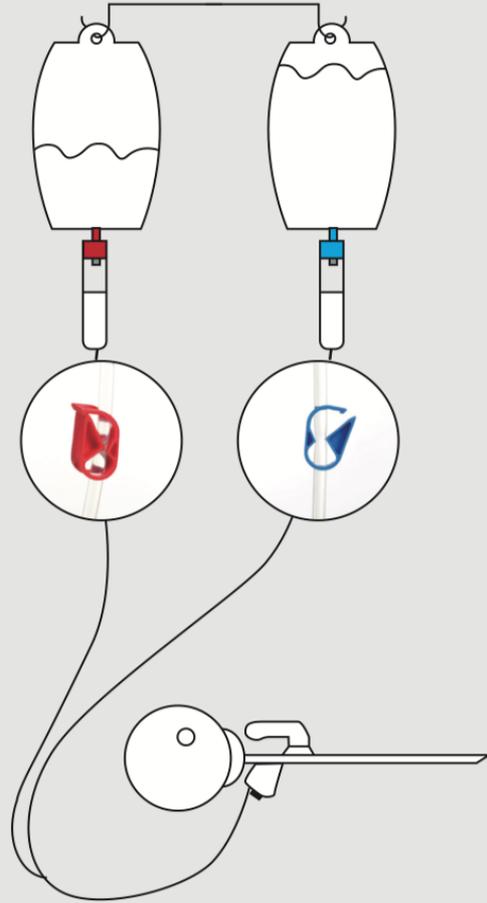


Double barrel tubing can be easily separated after the procedure and clamps conveniently moved just under the drip chamber to facilitate post-operative irrigation

REF	DESCRIPTION	QUANTITY
URO 005	NO BUBBLE TURP	24



Wismed CC  
71 Princess Alice Avenue,  
Durban, 4001, South Africa  
Tel: +27 31 206 0177  
orders@wismed.co.za  
www.wismed.co.za

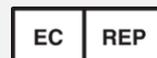

TNMC Devices (Llc)  
Ulica Franje Krezme 1  
10000 Zagreb, Croatia  
Tel: +38 59 1622 9658  
Email: tnmchr@tnmc.biz